



ANLAN

ULTRASONIC SKIN SCRUBBER

USER MANUAL

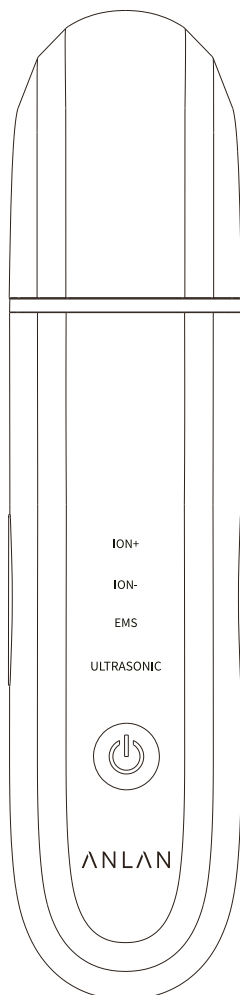

Overview

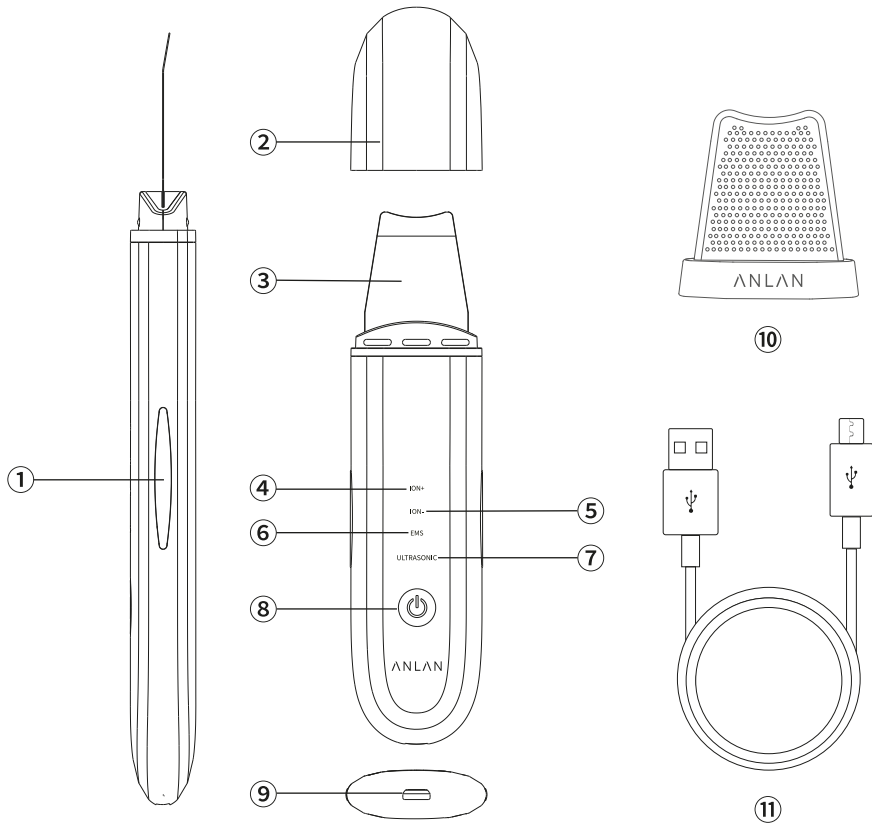

① Ion Conducting Bar

② Protective Cap

③ Spatula

④ Ion Export

⑤ Ion Import

⑥ EMS

⑦ Ultrasonic Care

⑧ ON/OFF & Mode Button

⑨ Charging Port

⑩ Silicone Cover


⑪ Type-C Cable

Specifications


Name: Ultrasonic Skin Scrubber
Materials: ABS
size: 171×40×14 mm
Weight: Appr. 75g
Charging Time: Appr. 2h
Battery Capacity: 500mAh
Standby: 50~60min
Rated Voltage: 3.7 -5 V
Automatic Shutdown Time: 5 minutes/mode
Accessories: Type-C Cable, Storage Bag, User Manual, Silicone Cover
Warranty: 1 year




Ultrasonic Skin
Scrubber*1

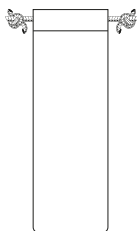


Storage
Bag*1

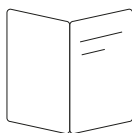

User
Manual*1

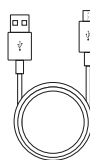


Charging
Cable*1

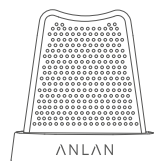


Silicone
Cover*1

Preparing Before Using

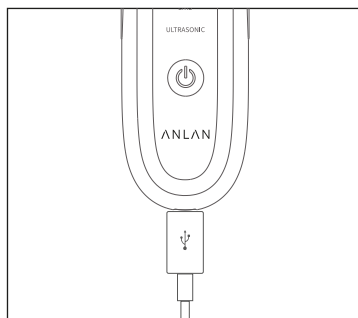

- ① Fully charge the scrubber, LED INDICATOR will keep on after fully charged.

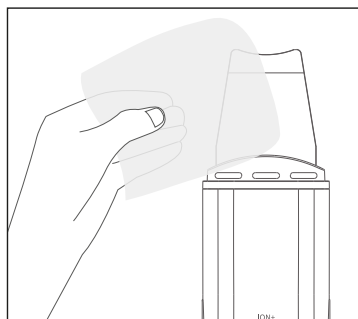


- ② Clean the spatula with wet tissue or wipes.

How to Use



- ※ Short press to turn on the machine and after power-on short press the power button again to change modes. Long press the power button to turn off machine.

Ultrasonic: Ultrasonic vibration

- ① Wash your face , then wipe your entire face with water or makeup remover.
- ※ Never use the device on dry skin.

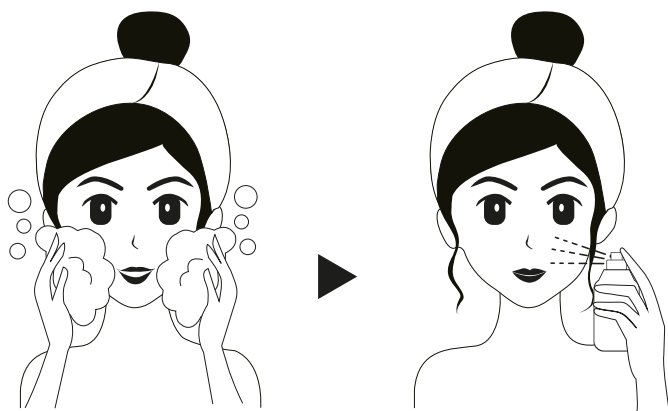


- ② Press the on/off button to start the machine. Briefly press the button and select the ultrasonic cleaning mode. Hold the ion conducting bar, carefully attach the spatula to the skin at an angle of about 30-45 ° and move it over the face in the direction of the arrow.
- ※ The spatula cannot be used on the same part of the skin for a long time and must be moved constantly when it is used.
 - ※ Depending on the brightness of the room, the water mist generated by the vibration of the spatula head will be different.

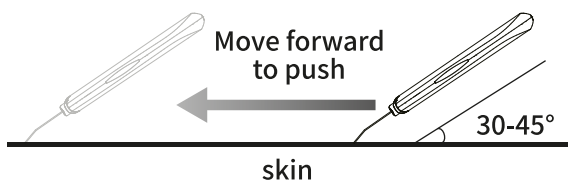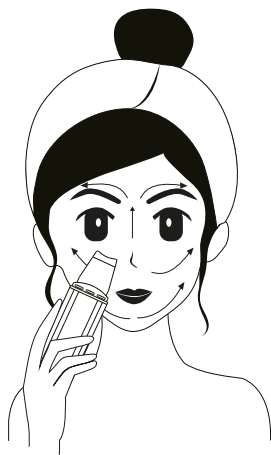


- ③ After use, wash your face to remove the sebum and oil, then wipe the spatula off with a towel or cloth.
- ※ Only the spatula head is waterproof, please do not wash other parts of machine.

Ion+: Ultrasonic vibration + Ion export

- ① After washing, wet the face skin with water or makeup remover.
- ※ Never use the device on dry skin.

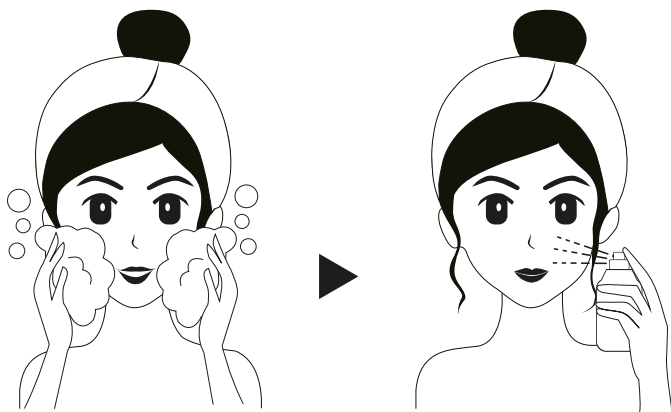


- ② Press the on / off button to start the machine. Briefly press the button and select the ion+ mode. Hold the ion conducting bar, carefully attach the spatula to the skin at an angle of about 30-45 °and move it over the face in the direction of the arrow.
- ※ The spatula cannot be used on the same part of the skin for a long time and must be moved constantly when it is used.

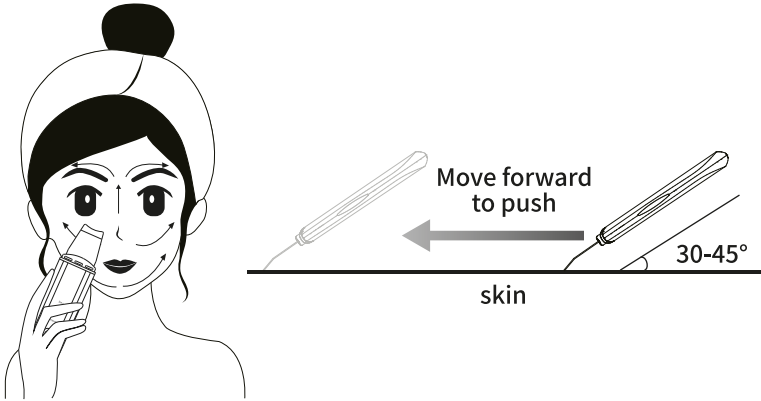


- ③ After use, wash your face and spatula to remove the sebum and oil, then wipe it with a towel or cloth.
- ※ Depending on the brightness of the room, the water mist generated by the vibration of the spatula head will be different.

Ion-: Ultrasonic vibration + Ion Import

- ① After removing makeup and washing your face, apply mask or moisturizer on the skin.
- ※ After using the ION - mode, it will facilitate the penetration of the beauty ingredients.

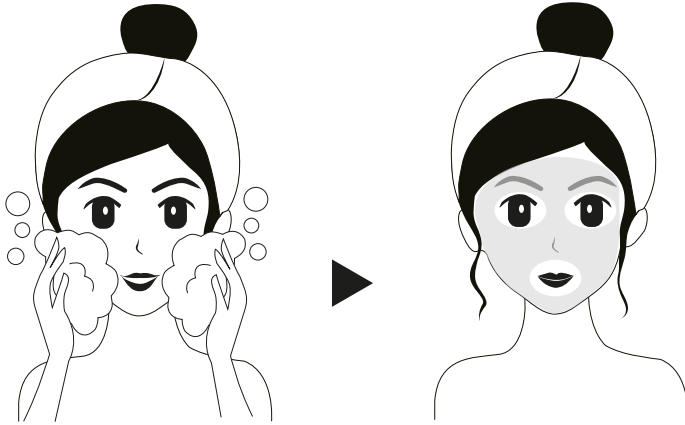


- ② Press the on / off button to start the machine. Briefly press the button and select the ION- mode. Hold the ion conducting bar, carefully attach the spatula to the skin at an angle of about 20-30 ° and move it over the face in the direction of the arrow.
- ※ The spatula cannot be used on the same part of the skin for a long time and must be moved constantly when it is used.

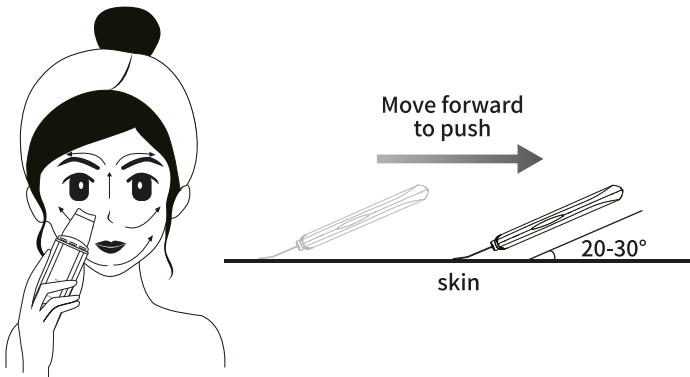



- ③ Use creams to nourish your skin. Wipe the creams on the head of the device with a paper towel.
- ※ Only the spatula head is waterproof, please do not wash other parts of machine.

EMS Lifting : Ultrasonic vibration + Ion Import + EMS

- ① After cleaning in export mode, wet the face skin with serum, essence or cream.
- ※ Please never use it on dry skin.

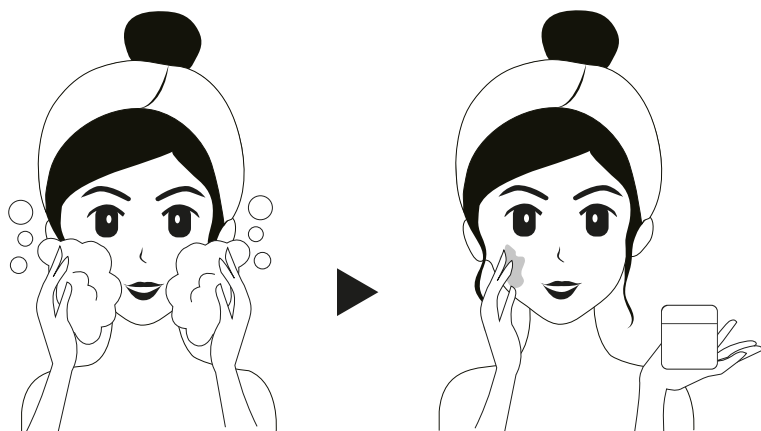




- ② Press the on/off button to start the machine. Briefly press the button and select the EMS mode. Hold the ion conducting bar for clear current, carefully attach the spatula to the skin at an angle of about 20-30 ° and move it over the face in the direction of the arrow.
- ※ The spatula cannot be used on the same part of the skin for a long time and must be moved constantly when it is used.

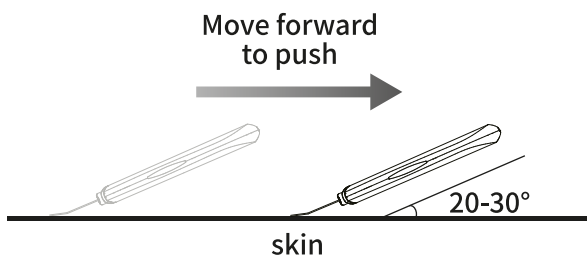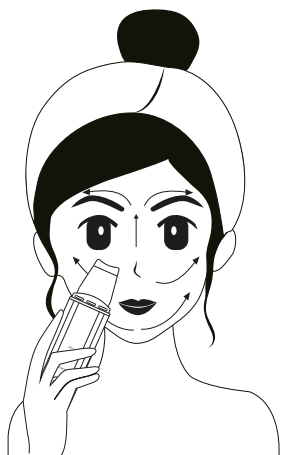


- ③ After use, use the beauty serum to condition the skin. Wipe the creams on the head of the device with a paper towel.
- ※ Parts other than the spatula head cannot be washed by water.

Attentions while using:


- Follow the directions to move the Spatula.

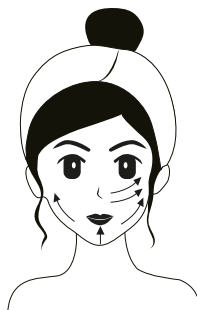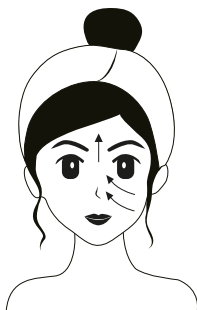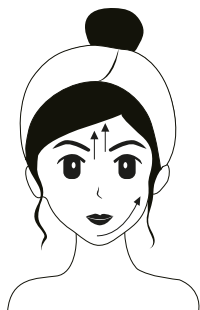

- Attention to the direction of the spatula while using:

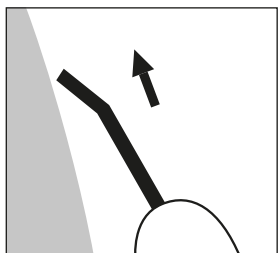


Pic.1:
Ultrasonic/ Ion + Modes

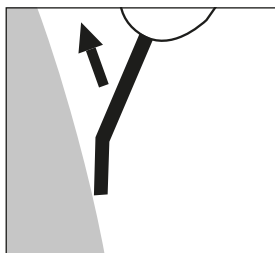


Pic.2:
Ion- / EMS Mode


How To Clean Device

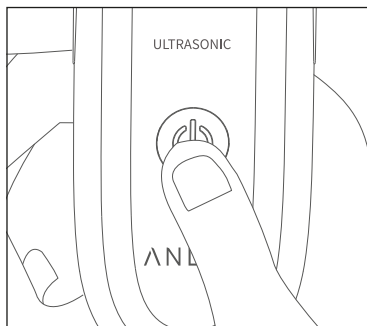

- ① After using, long press the power button for 3-5 seconds to turn off the device.

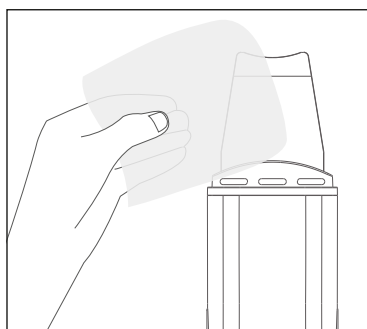

- ② Wipe off the liquid or essence remaining on the device.

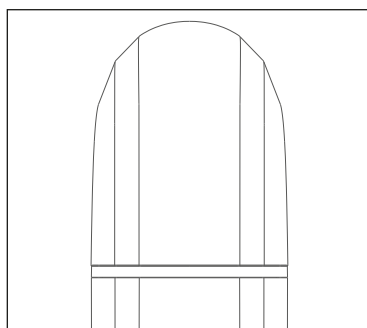

- ③ Clean the spatula with water or a cotton pad, and cover the protective cap after it is dry.

Precautions

- Remove your makeup before using, and keep your skin moist enough during use.
- Keep the SPATULA clean before and after use. Clean it with wet cloth or cotton pad. Do Not wash it with water directly.
- Do not soak the device in water, or use it while bathing.
- Do not use the device on the same area of your face or neck for a long time.
- Do not use it on eyelid or other parts but face.
- Do not use it for any purpose other than as described.
- You are recommended to use the device 1~2 times a day, no more than 15 minutes each time. The same part of your face or neck should be used again after 6 hours later.
- Pregnant women or women under menstruation period are prohibited.
- Heart disease, high blood pressure, coronary heart disease, current allergies are prohibited.
- Infectious skin diseases, burn skin diseases, anticoagulants, or skin allergies are prohibited.

Q: It may be more effective to press firmly?

A: No, please don't press hard to avoid straining the skin.

Q: What kind of cream should I use?

A: You can use it with regular creams (essences/lotion/high hydration facial creams are fine).

Q: Is it effective for men too?

A: This product is unisex. Each mode looks great on men's skin.

Q: Is it okay to use every day?

A: Ultrasonic and Ion + can be used 2-3 times/ week, EMS and Ion- can be used everyday.

Guarantee

First of all, we hope you never need this warranty.

But we also want to offer you the best service, especially if the product has problems.

Our goal is to offer our customers the best products with first class quality and service. We guarantee return and modification within 30 days, and technical support at any time.

If you are not satisfied or have any questions about the ANLAN product, please do not hesitate to contact us.

service@anlan.co

Thank you,

ANLAN customer service

Q&A

Q: Is it waterproof? Can I use it in the bathroom?

A: Just the shovel head is waterproof, the body is not waterproof, don't rinse it with water. you only need to clean the shovel head with a clean paper towel. Do not use or place the product where it is exposed to water (such bathroom or damp place).

Q: The device automatically shuts down suddenly, and the body became very hot.

A: Check if there is no power, you can turn off the power and charge it. The device has a 5-minute automatic shutdown function. Ultrasonic vibration can cause heating of the head and body ends, which is normal.

Q: The hands will feel tingling or pain during use, Is it a leakage of electricity?

A: This is not a leakage, the product through hold the Ion conducting bar with fingers to form a current loop, which will produce EMS micro-current to stimulate muscle massage, reduce atrophy, remove wrinkles, anti-aging, lift tighten your face, improve skin lines to restore skin elasticity. Some people may feel uncomfortable. Please be assured that micro currents are harmless to the human body. If you do not like this feeling, please do not touch the Ion conducting bar.

Q: Can blackheads be removed?

A: High frequency ultrasonic vibration can remove dirt from pores to reduce acne, keratin and blackheads. (Note: If you encounter stubborn blackheads, acne, whiteheads, we recommended to use it with blackhead remover tool)

Troubleshooting

If the device is not working after pressing power button:

- No battery power. Restart again after fully charged.
- Not working after fully charged. Contact customer service for repair or replacement.

If the device cannot charging:

- Charging cable is not well connected or the charging cable is damaged. Check the charging connection or use another cable.
- If it still does not work. Contact customer service for repair or replacement.

If the machine beeps three times:

- Low battery. Restart again after full charged.
- If it's not working after fully charged. Contact customer service for repair or replacement.