


ANLAN

ANLAN Body Beauty Device User Manual

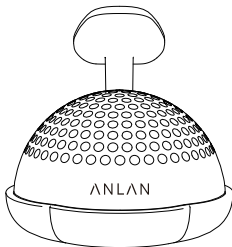

※Thank you for purchasing ANLAN product. This product is not a medical device used for diagnosis and treatment, but a support product for general household use.

Be sure to read the instruction manual carefully before use to ensure correct and safe use. In addition, this document also serves as a warranty card, so please keep it in a safe place.

Safety Precautions

In order to prevent possible injury and property damage to users from improper use, please read this manual carefully and follow the manual instruction.

 Danger	Inappropriate use can lead to injury or property damage.
 Warning	Mishandling may result in death or serious injury.

Warning

Do not use the following medical electrical equipment with this device together:

- Implantable medical electronic devices such as pacemakers
- Medical electronic devices for life-threatening items such as heart-lung machines
- Wearable medical electronic devices such as Electrocardiographs

The following people do not use

- People with heart disease
- People with skin allergies
- People who cannot sense temperature

- People with facial neuropathy disease
- People with skin diseases and atopic dermatitis
- People with blood abnormalities such as bleeding disorders and hemophilia
- Toddlers, growing children, and people who cannot express their intentions
- People with sequelae on their heads due to traffic accidents, etc.
- People who are pregnant and within one month after giving birth

Do not use in the following cases or body parts

- If you are not feeling well
- Lips, eyelids, ears, neck area with arteries
- Skin areas that have been depilated on the day
- Skin areas with scratches, eczema, swelling, acne, etc.
- Sites undergoing surgery (plastic surgery, cosmetic surgery etc.)
- If you have an abnormality on your skin such as dermatitis or excessive sunburn
- When taking steroid hormones, immunosuppressive drugs, or drugs that cause drowsiness


The following people, please be sure to consult with your doctor before using.

- People with acute illness
- People with febrile illness
- People with infectious diseases
- People with malignant tumors
- People with extremely sensitive skin
- People with a history of serious illness
- People with sensory impairment due to severe peripheral circulatory disorders such as diabetes
- Patients who need postoperative rest or those who have experienced physical discomfort due to surgery within two weeks

- Do not put device into high temperature places such as fire,
- Do not modify or repair it to avoid possible fire, electric shock, or injuries
- Do not disassemble by yourself to avoid possible fire, electric shock, and injury
- Keep it in places out of the reach of children to avoid possible accidents or injuries .
- Do not charge in humid place like bathroom to avoid short circuit and electric shock
Do not insert or unplug the adapter from socket with wet hands to avoid electric shock, or injury

- Do not charge, use or place it at high temperature to avoid built-in battery from heating, inflaming, or exploding

Product Details

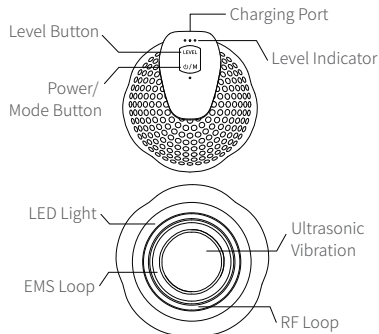

Pack List



Main Device* 1



Charging Cable* 1



User Manual* 1



Silicone Base*1



Storage Pouch*1

Specifications

Product Name	Body Beauty Device
SKU Model	OFY-2102
Size	Appr. 85*85*84mm
Weight	Appr.170g
Material	ABS
RF	1.9MHz
EMS	80KHz

Ultrasonic Wave	3.2MHz
Light Wavelength	Red 620±10nm Blue 470±10nm Purple 390±10nm
Vibration Motor RPM	Appr. 9000 RPM
Charging Port	Micro USB Charging
Battery Capacity	1100mAh
Working Time	Appr.40 minutes (Max Power)
Charging Time	Appr. 180 minutes
Input	DC 5V $\overline{\text{---}}$ 1A
Voltage	7.4V
Smart Timer	10 minutes
Waterproof	IPX6
Warranty	1 Year


※The actual working time will vary depending on the conditions of use and operating environment.

Before Use

Please do a skin patch test:

Test with the skin care products you usually use to make sure its compatibility with the machine to prevent possible allergies.

※Please be sure to do a skin test to prevent any problems.

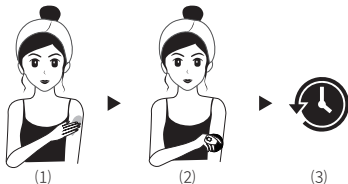

1. Wash your arm with soap and wipe it with a towel.
2. Turn on the product, select mode 1, do skin test by using it with your daily skin care products
3. If there is no abnormality after 24 hours, you can rest assured to use it.

※If any abnormalities (such as intense itching or pain) occur after patch testing, the product may not be compatible with the product. Please stop using the product.

How To Hold The Device

Hold Device:

Hold device the way above

How to apply to the skin

Make the holding head perpendicular to your skin and put entire flat metal surface in contact with your skin. Tilting the head may cause incomplete contact and skin abnormality.



About Buttons

⏻/M Power/Mode button

Press and hold button for 3 seconds to turn the power on and off.

Press 1 second to switch between the three modes.

Mode 1 → Mode 2 → Mode 3 (repeat).

LEVEL Level button

Short press to change between three levels or 3 vibration mode. Level 1→level 2→level 3 or vibration 1→vibration 2→vibration 3. (Repeat)

About All Modes

Mode	Functions Combo	Use Frequency	Suggest Time
Mode 1	Ultrasonic+ Vibration+Blue Light	2-3 times/ week	10min- utes
Mode 2	RF+EMS +Red Light	2-3 times/ week	10min- utes
Mode 3	Ultrasonic+ Vibra- tion+RF+EMS+Purple Light	2-3 times/ week	10min- utes

※You don't have to use all modes in one day. Please choose proper mode and control using time according to your skin condition.

※Caution: Overuse can stress your skin, so avoid using all modes on the same area every day.

※After 10 minutes of power on, the power will be turned off automatically. Switching modes midway auto-off time will count down again.

Usage Notice

※It is recommended to use after bathing.

※We recommend using it with water-soluble lotion or gel.

※Be sure to use by moving on the skin. Stay on the same area too long may result in burns.

※Frequency of use: depends on your own need and feeling, such as one area for 10 minutes or 5 minutes on each of 3 areas. 3 times a week is recommended.

※After 30 minutes, please moisturize your skin with lotion.


How To Use:

※Do not use on face.


1. After bathing as usual, apply lotion or serum to the skin.
2. Long press the power/mode button to turn on the machine, and short press to select Mode 1/Mode 2/ Mode 3.
3. Short press the level button to adjust the level. Adjust the level from 1 to 3 according to your skin condition.

- ※Be sure to try from level 1 on first use and change it until you get used to it.
4. Hold the main body with dry bare hands, and use it while sliding, so that the entire electrode head touches the skin and does not stay in one place for too long.
 5. Press and hold the Power/Mode button to power off.


Recommended Usage For Body

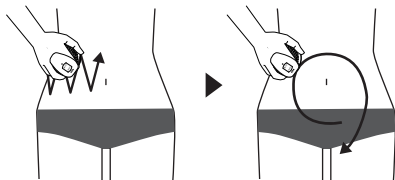

- Abdomen:Moves device up and down in a zigzag on belly and make circular motions around your navel.

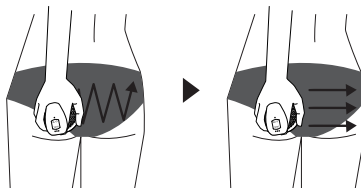

- Buttocks: Moves device up and down in a zigzag and push from the inside to the outside of the buttock

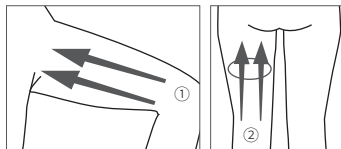


- Thighs:
 1. Massage from the inner knee to the inner thigh near the groin
 2. Push device from the back of the knee upwards to the buttocks to massage



- Calf: Move from back of ankle to back of knee to massage

- Arms: With palms facing up, lower arms straight and massage from inner elbows to close to armpits.

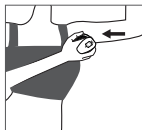

About Charging

- During Charging: When the adapter is connected, the indicator light flashes.
- Full charge: Approximate charging time is about 3 hours. When fully charged, the light stays on.
- Low battery: the light flashes all the time until it's dead.

- When using for the first time or not using for a long time, please charge it before use. In addition, in order to protect the built-in battery, please charge it once a month even if it is not used for a long time.
- The main unit may become slightly warm while charging, but this is not a malfunction.
- To prevent overcharging, please unplug the cable after charging is complete. If you continue to charge the battery for a long time after full charging, the built-in battery will be exhausted and its life will be shortened.
- This product cannot be used while charging.
- Please make sure the charging port is clean and dry before charging.
- Be sure to use the voltage and AC adapter specified in this document.

Maintenance

- Turn off the power after use in time.
- The main body cannot be washed with water.

- The surface of the head is IPX6 waterproof, but the charging port and buttons are not waterproof, so please do not put it in water or splash it with water. If there is dirt, wipe it off with a wrung wet towel or paper towel.
- Never use thinner or benzine.
- If it is dirty, wipe it off with a towel or tissue.
- Do not put the main body in a fire or overheat.
- Do not leave product in direct sunlight or near heat .
- No one but a repair technician should disassemble, repair or modify.
- Do not leave the product in a humid place such as bathroom.
- Store it out of reach of children and in a safe place.
- Please do follow precautions and warnings above.
- We are not responsible for any injuries, malfunctions, damages caused by your incorrect use, disassembly, modification, repair, etc.

Trouble Shooting

Problems	Possible Causes	Solutions
can not turn on	You haven't long press power button for seconds.	Long press power button for seconds
Only work for a few minutes after charging	Insufficient charge	Please fully charge it
	Product failure	Please contact ANLAN Customer Service for help
The device is hot	Used continuously for too long	Turn off power and rest for 10 minutes or more.
Doesn't get cold in COOL mode	Used at improper temperature	Please use in an environment of 15 ° C to 30 ° C.

Problems	Possible Causes	Solutions
Ion / EMS is strong (painful)	Set improper level	Choose lower level
	Device top not touch skin wholly	Make sure whole top touch skin
Ion / EMS is weak No feeling of EMS	Set improper level	Choose higher level
	Your skin is too dry to transfer EMS Not perfectly touch skin	Moisten your skin with lotion before using and make sure to touch skin entirely
RF feels too hot	Device top stay too long on same part of skin	Move device to massage skin
	Move device too slowly	Move device faster

Problems	Possible Causes	Solutions
No feeling of RF	Your skin is too dry	Moisten your skin with gel before using
	Move device too fast	Move device more slowly

FAQ

Q: I feel a little tingling on my skin when I use it, is it normal?

A: When using mode with EMS, it is normal to have a slight tingling sensation.

You can will feel more comfortable to do this

①Keep skin moisturizing enough with special gel.

②Use it from level 1 at the beginning, and gradually increase it as you get used to it.

Q: How often can I use it?

A: It's recommended to use 10 minutes in the same area once and 3 times a week

※ The above is for reference only. Please use it according to your skin condition and preference.

Q: Can it be used on the face?

A: Better not. This device with higher energy than ANLAN other face massager and may put too much pressure on sensitive face with low fat. We recommend using this product on areas such as arms, calves, thighs, belly, hips and buttocks.

About Disposal of This Product



This product uses a lithium-ion battery. Lithium-ion batteries are a valuable recyclable resource. When disposing of used products, remove the lithium-ion battery and take it to a cooperating store that has a recycling box or a specific collection window.

Warranty & Guarantee

First of all, we hope you never need this warranty.

But we also want to offer you the best service, especially if the product has problems.

Our goal is to offer our customers the best products with first class quality and service. We guarantee 365 days warranty for ANLAN products.

If you are not satisfied or have any questions about the ANLAN product, please do not hesitate to contact us.

service@anlan.co