


ANLAN

Facial Beauty Device

User Manual

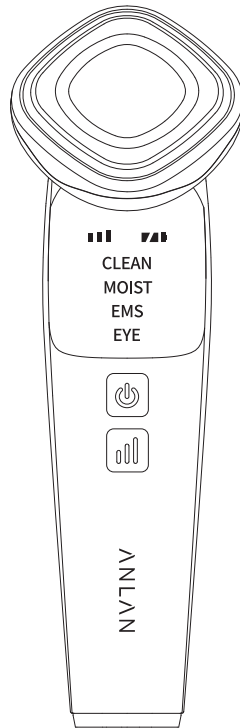

Thank you for purchasing the ANLAN Facial Square Beauty Device. Please be sure to read this manual before use to ensure safe and correct use. In addition, this manual also serves as a warranty card, so please keep it in a safe place.

※This product is not a medical device used for diagnosis or treatment, but a support product for general household use.


Safety Precautions

In order to prevent harm to users, please read this manual carefully and follow the manual step.

 Warning	Inappropriate use can lead to undesirable consequences.
 Attention	Please follow the instructions to use the machine to avoid hurting.

 Prohibition	 Warning
	• Never use this device with following medical electrical equipment. • Tiny electro-magnetically damaged implant, such as pacemaker. • Medical electronics for life support such as artificial heart / artificial lung. • Portable medical electronics, such as device of electrocardiogram.
 Mandatory	 Attention
	The people mentioned below should consult a doctor before using the device. • Women in pregnancy • People with acute (painful) illness • Patients with tuberculosis • Patients with facial neuralgia • People with heart problems • People with malignant tumors • People with abnormal blood pressure • People suffering from infectious diseases • People with hemophilia • People with febrile illness • People with allergies • People suffering from infectious diseases • People with a history of dermatitis like atopic and contact dermatitis • People with sensory disturbances due to severe peripheral circulatory disturbances such as diabetes

Product Details

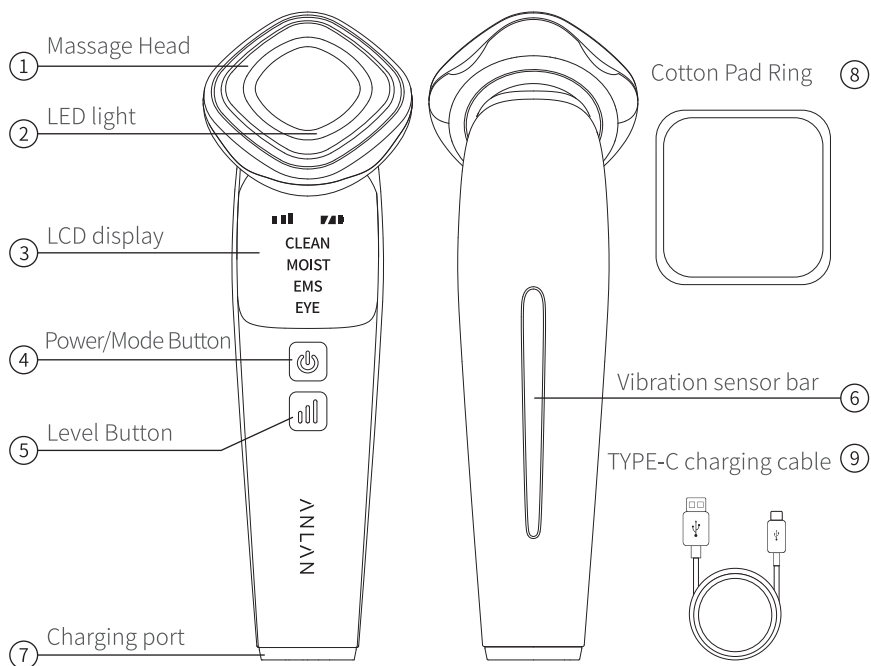


Specification

Product Name	ANLAN Facial Beauty Device
Model	D835
Material	ABS
Product Size	170*55*56mm
Product Weight	About 115g
Vibration Frequency	About 10000rpm
Temperature	Hot: About 38 °C ~ 42 °C
Wavelength of Light	Red Light About 623nm, Blue Light About 463nm
Battery Capacity	450mAh
Charging Time	About 2.5 Hours
Automatic Power OFF	About 5 minutes/time, switch mode and re-time
Accessories	TYPE-C charging cable, Pad Ring, user manual

Before Use



For the first use, we recommend you test device on your arm. Try it on the skin with your toner, lotion or cream before long term using.

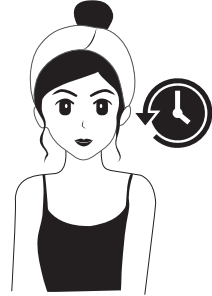

1. Wash your arms with facial cleanser, then wipe with a towel.

2. Select the "Firming and Lifting" mode and test it with your skin care products.

3. If there is no discomfort after 24 hours, you can use it with confidence.



Note

Remove makeup and cleanse face before use.

- Do not stay in one place for more than 10 seconds.
- Please adjust a suitable temperature according to your needs.
- Be sure to touch the vibration sensor bar when using it.
- Do not use directly on dry skin, it is recommended to use skin care products.
- If you feel unwell, please stop using it and seek medical advice in time.

Modes & Functions

Mode	Function	Recommend	Use Frequency	Level
CLEAN	Heating + red light + positive ion export + access to Skin vibration	Use with cleansing water,cleansing cream, toner, etc.	Once a day, 5 minutes / time	Temp:3 Level (38°C-40°C-42°C)
MOIST	Heating + blue light + negative ion introduction + access to Skin vibration	Use with mask and skin care products	Once a day, 5 minutes / time	Temp:3 Level (38°C-40°C-42°C)
EMS	EMS+red and blue light alternating	use with gel	Once a day, 5 minutes / time	EMS:3 level.
EYE	EMS hot compress Blue/ red light	Use with lotion, milky lotion, cream etc.	Once a day, 5 minutes / time	Temp: 3 Level (38°C-40°C-42°C) EMS: 3 Level


※You don' t have to use each mode in order. You can choose the mode and usage time according to your skin condition.




Note

Excessive use will put a strain on your skin, so please avoid using all modes every day.


How to Attach / Detach the Cotton Pad Ring

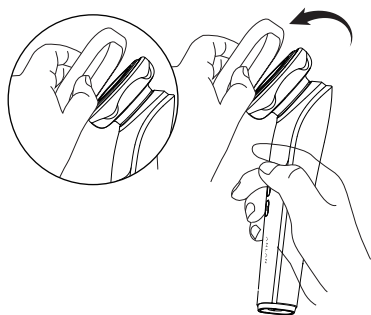


- 1** Remove the cotton pad ring.

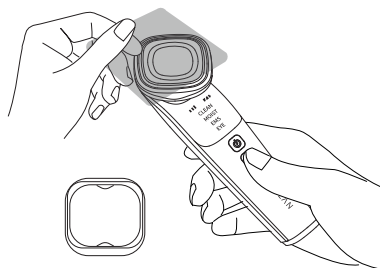

- 2** Put cotton on the head.

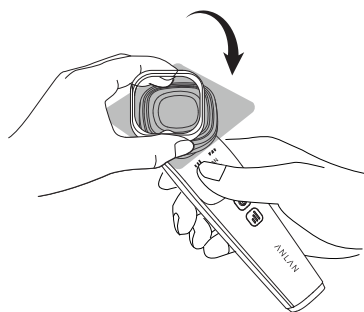

- 3** Install the cotton pad ring.

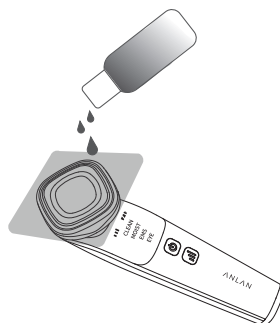

- 4** Apply makeup remover/toner to cotton.

How to Use



Power / mode button


Press and hold 2 minutes to turn on / off.
Short press to adjust the modes.



Level button

Short press to adjust temperature and EMS intensity.

※ When using it for the first time, please be sure to use it from level 1 until you get used to it.

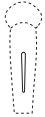


Vibration sensor bar

When using the vibration function, the hand must touch the vibration sensor bar on the back of the fuselage.

Using Steps


1. The deep cleansing mode works with makeup remover, makeup remover oil, makeup remover cream and other products to remove makeup.

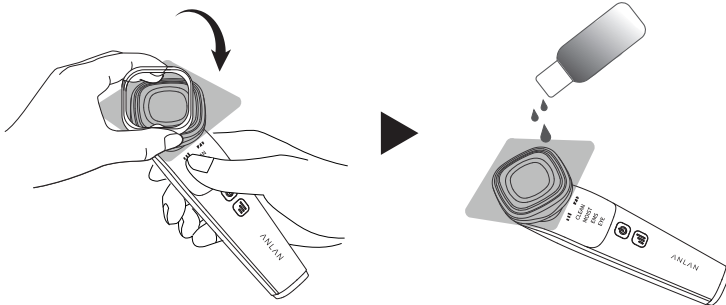


2. Wash your face with facial cleanser and wipe off the moisture.



3. The deep cleansing mode is paired with toner for secondary cleansing.

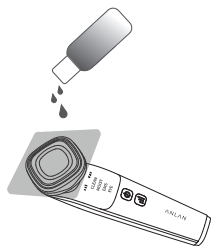

4. When using nourishing introduction, firming and lifting, and eye care, it can be used with lotion, essence, gel, mask, etc.

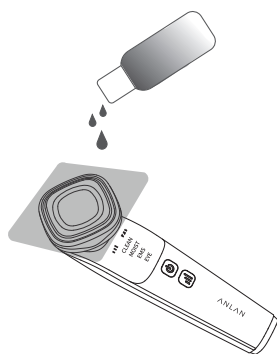

5. After using skin care products, move the massage device slowly.

Hold the vibration sensor bar of the device, if the hand does not touch it, the product will not vibrate.

- ① From nose to forehead
- ② From the bottom of the forehead to the top
- ③ From the inner side of the face to the outer side
- ④ From chin to cheek



Maintenance

- Turn off the power after use.
- The main body cannot be washed with water.
- Never use thinner or benzine.
- If it is dirty, wipe it off with a towel or tissue.
- Do not put the main body in a fire or overheat.
- Do not leave product in direct sunlight or near heat .
- No one but a repair technician should disassemble, repair or modify.
- Do not leave the product in in a humid place such as bathroom.
- Store it out of reach of children and in a place.
- Please do follow precautions and warnings above.
- we are not responsible for any injuries, malfunctions, damages caused by your incorrect use, disassembly, modification, repair, etc.

About Charging

- Never charge the product in wet places, such as when the terminal is wet or in the bathroom, there is a risk of fire, electric shock and burns.
- Before charging, make sure the charging port is clean and dry.
- If the battery power is low, the battery icon will blink and then shut down.
- Estimated charging time is about 2.5 hours. When charging the main unit, the battery level mark will blink. When fully charged, the battery The battery level mark will light up. (Note: Charging time may vary depending on charging environment.)
- Do not use while charging. Please unplug the cable from the main unit except when charging.
- Make sure to use the operating voltage, AC adapter, etc. under the conditions specified in this manual.
- If you are using this product for the first time, or you have not used this product for a long time, please charge it before use.
- Also, to protect the built-in battery, charge it once a month even when not in use.
- The room temperature suitable for charging is 15-30°C. Charging at low temperatures will shorten battery life, so please charge at an appropriate temperature.
- Also, charging speeds at low temperatures may drop. During the charging process, the main unit may become a little hot, which is normal and don't worry.
- If the operating time is unusually short even when the battery is properly charged, please contact customer service.

FAQ

Q: Does the head not vibrate?

A: ①To feel the vibration, the product must have a complete electronic circuit. Hold the vibration-sensing strip on the back of the device, and gently attach the massage head of the device to the skin to form an electronic circuit. (If the head is separated from the skin, it will not vibrate,because electronic circuits cannot be formed).

② Keep the skin moisturized with body lotion or essence. If the skin is too dry, the transmission of vibrations and EMS will be affected.

Q: Can the snap ring be removed and used?

A: The snap ring is only required for deep cleaning mode. The snap ring can be removed in other modes.

Q: Can it be used with makeup?

A: Please do not use it with makeup, please cleanse your face before use.

Q: My hands feel tingling during use, is it normal?

A: Please dry your hands before use. Since it works with a weak current, the current will also flow in the hand, forming a circuit. But this is a very small current, so it has no side effects on health.

Q: My skin has acne and inflammation, can I use it?

A: Since the skin is fragile when inflamed or acne-prone, please avoid sensitive areas.

Q: How long will it take to see the effect?

A: After 28 days of continuous use, you can see the obvious contrast effect before and after. It is recommended that users continue to use it as required, and take photos to record their transformation.

About the Disposal of Multifunctional Skin Care Machine



This product uses a lithium-ion battery. Lithium-ion batteries are a valuable recyclable resource. When disposing of used products, remove the lithium-ion battery and take it to a cooperating store that has a recycling box or a specific collection window.